

Software must be V3.0B or higher.


1. Install Tera Term application in your Laptop.

[See the sheet "How to install TeraTerm"](#)

2. Install a driver of Converter cable (Serial ⇔ USB)

Follow the setup procedure of Converter cable that you use.

3. Confirm that Converter cable is successfully recognized from Device manager

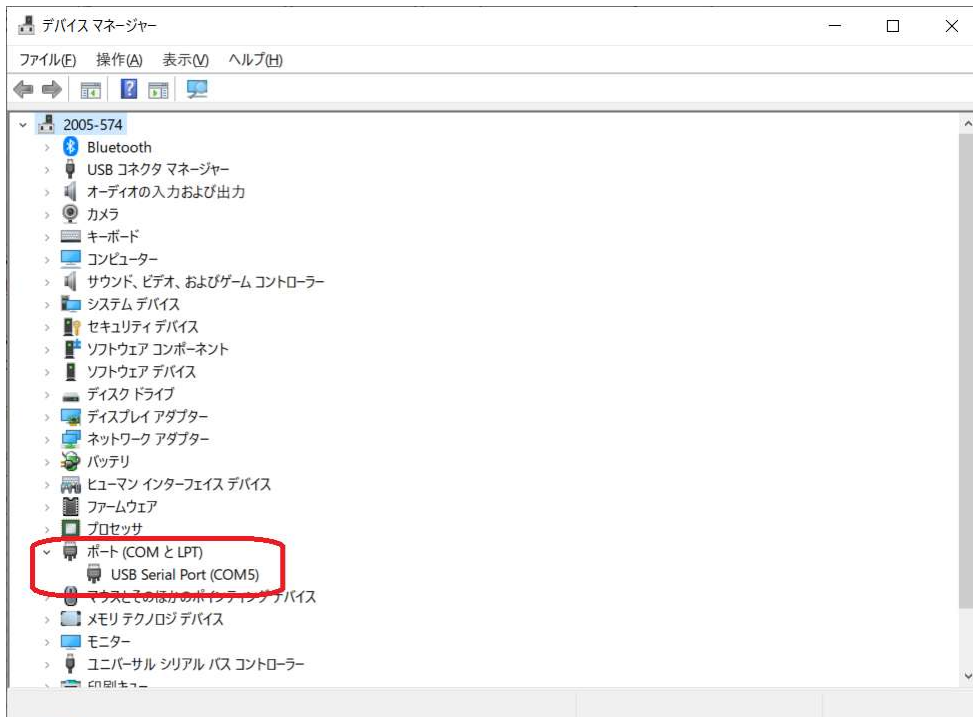


4. Go to adjustment mode of DBB-27



5. Change the first letter of serial number from "J" to "2"

Example: J6946405 → **2**6946405



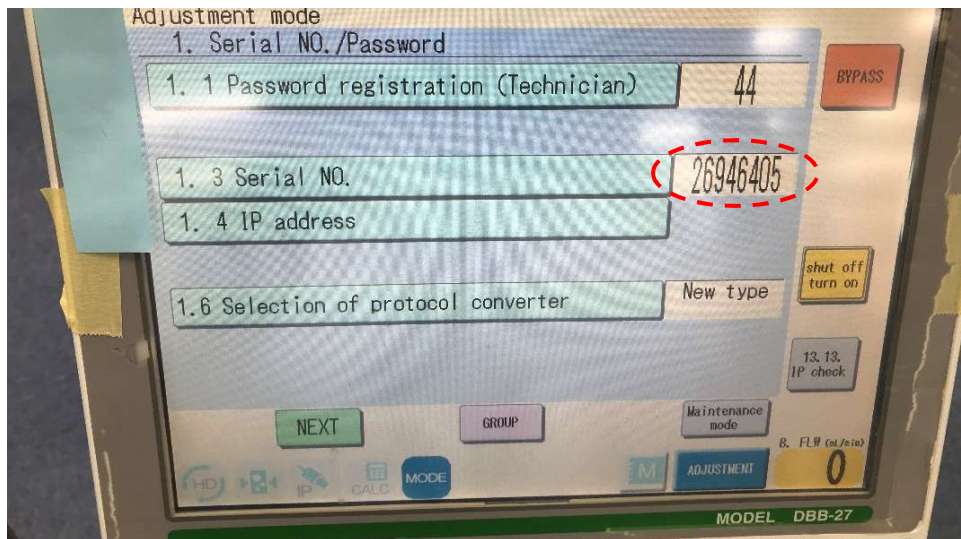

6. Change IP address to "0.0.0.0" and press "SAVE"

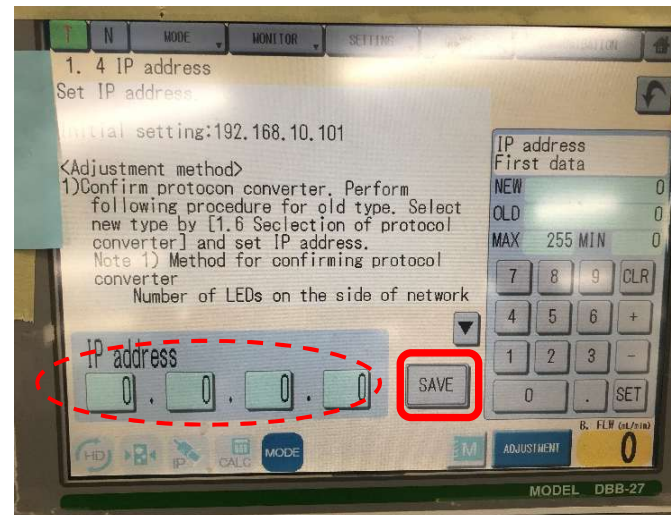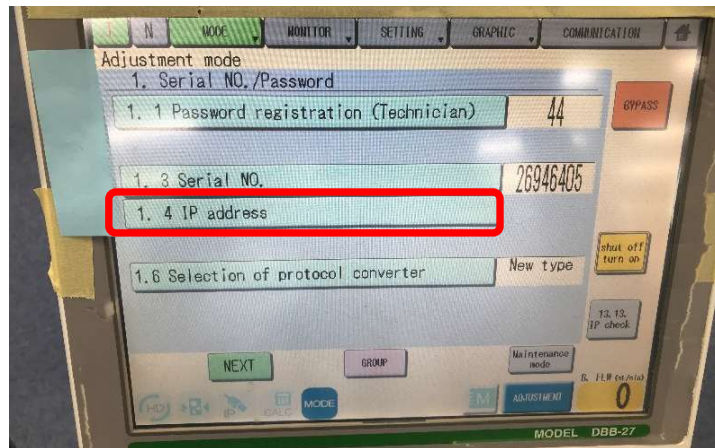


7) Set "Model selection of separate BPM" of 37. blood pressure monitor from Setting 2 to "0: built-in".

* If the hospital setting value is not "0: built-in" (use an external blood pressure monitor), record the setting value.

* To check the operation of external output 2, it is necessary to set "0: built-in".



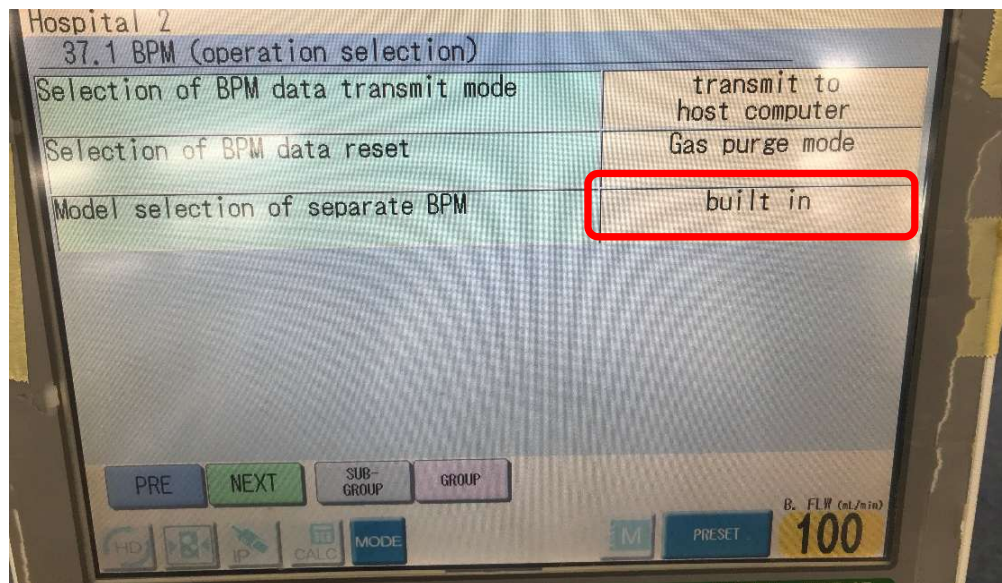


8. Turn OFF and Turn ON the machine.

IP address is stored to the machine after turned OFF.

9. Exit the adjustment mode.




10. Connect cable between PC and RS232C port2 of DBB-27

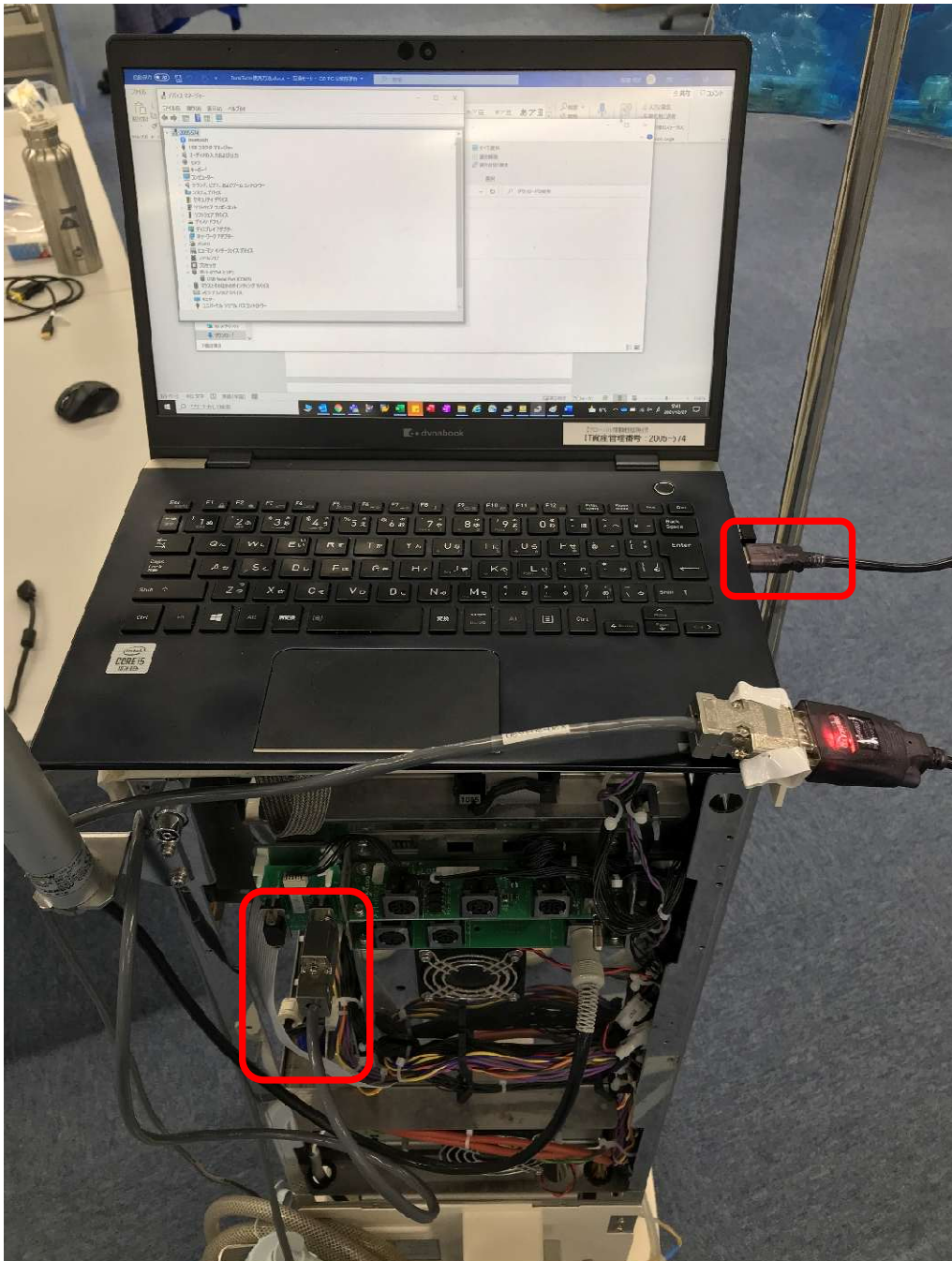



11. Open "TeraTerm" on your PC.



12. Select serial and confirm that correct port are selected.

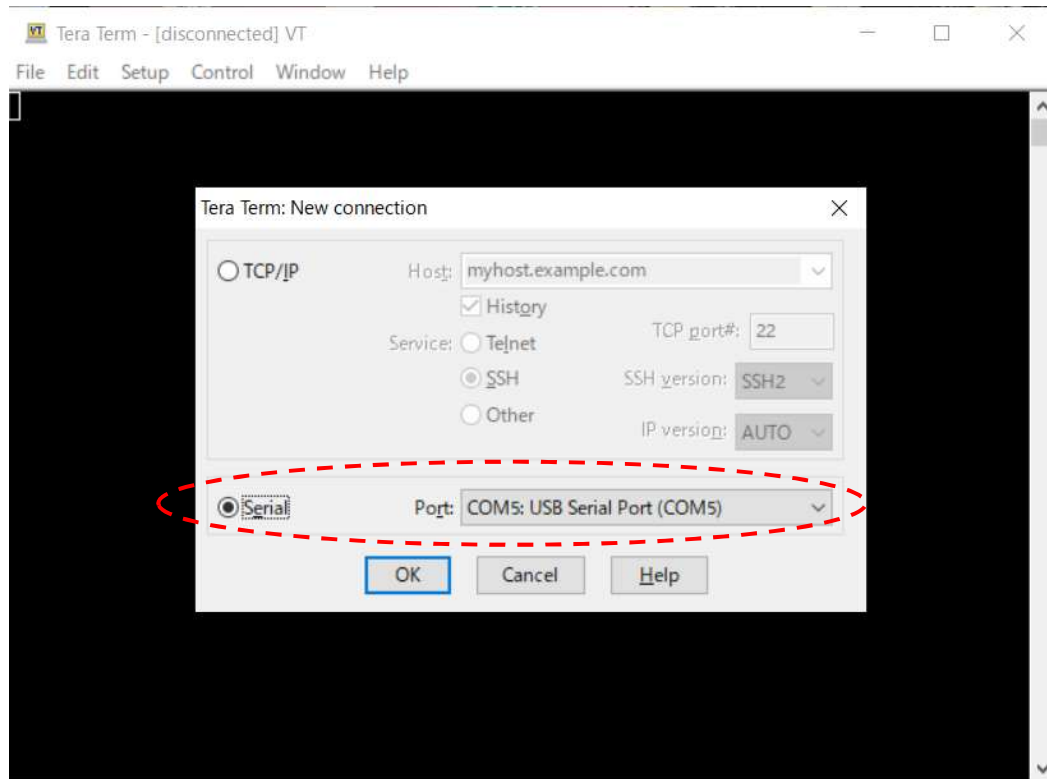



13. open "Serial port" window from "Setup" tab.

Change communication settings same as a picture below.

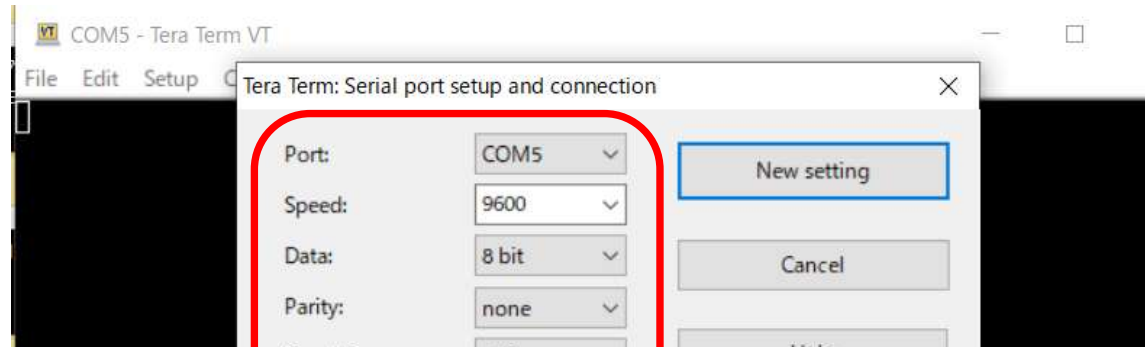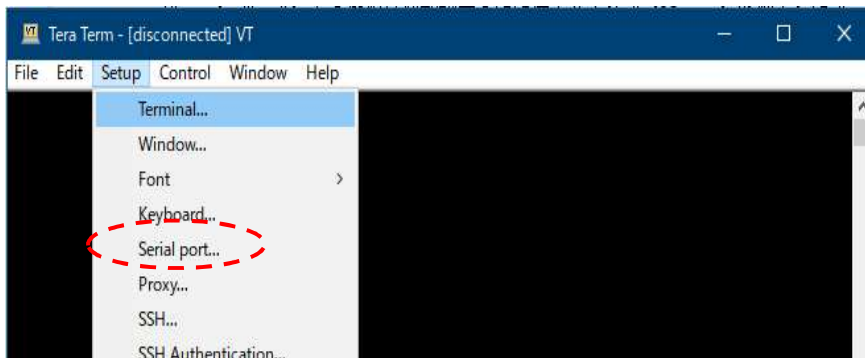

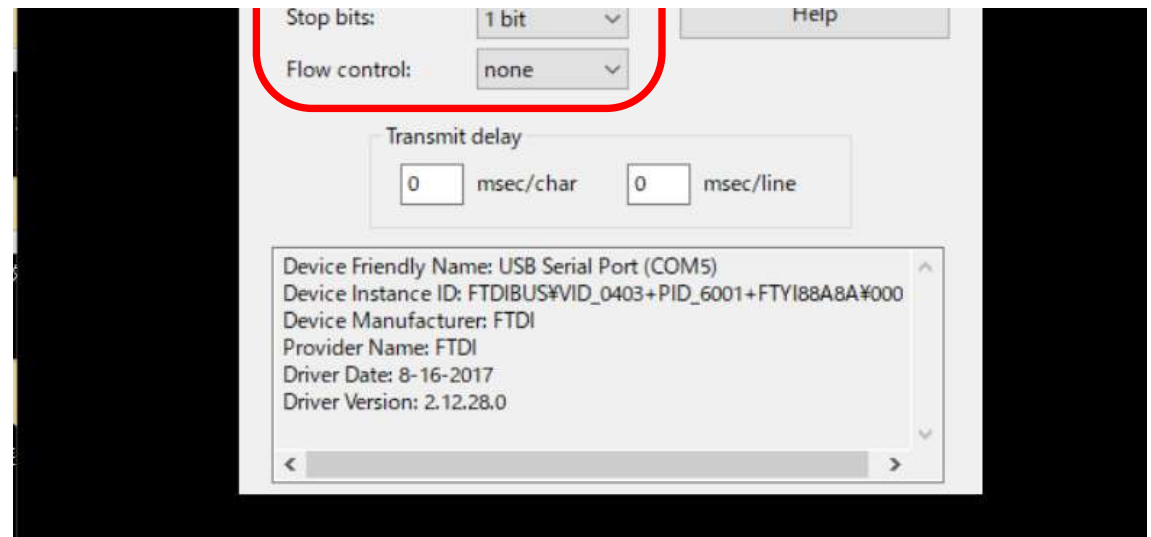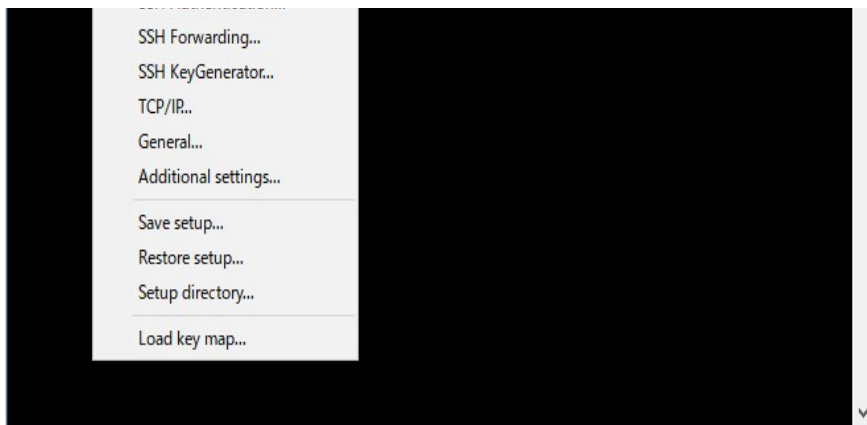


13. open "Terminal" window from "Setup" tab.
Change communication settings same as a picture below.

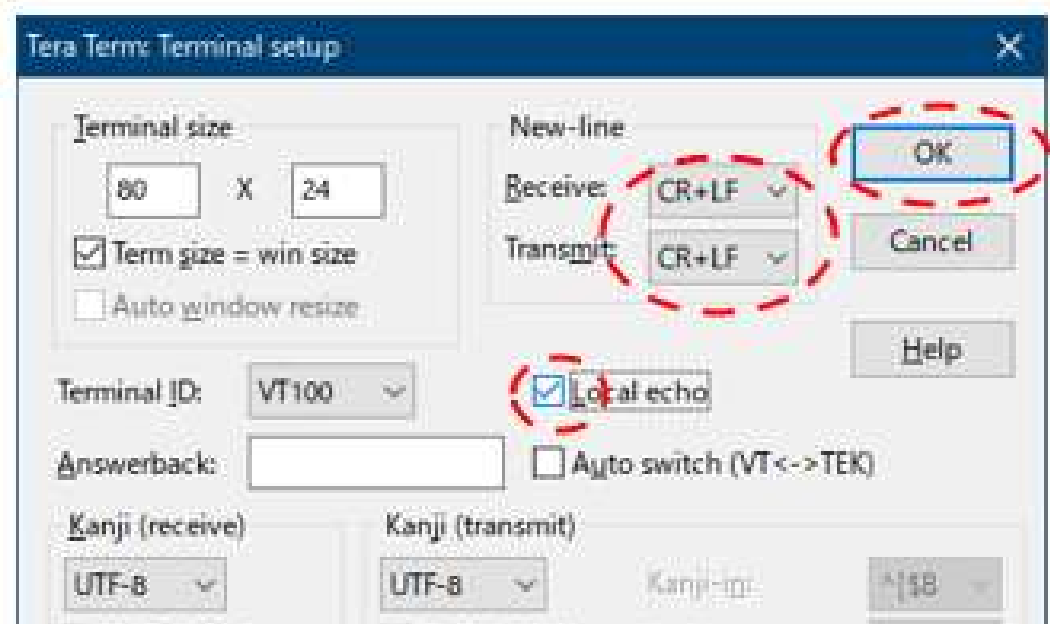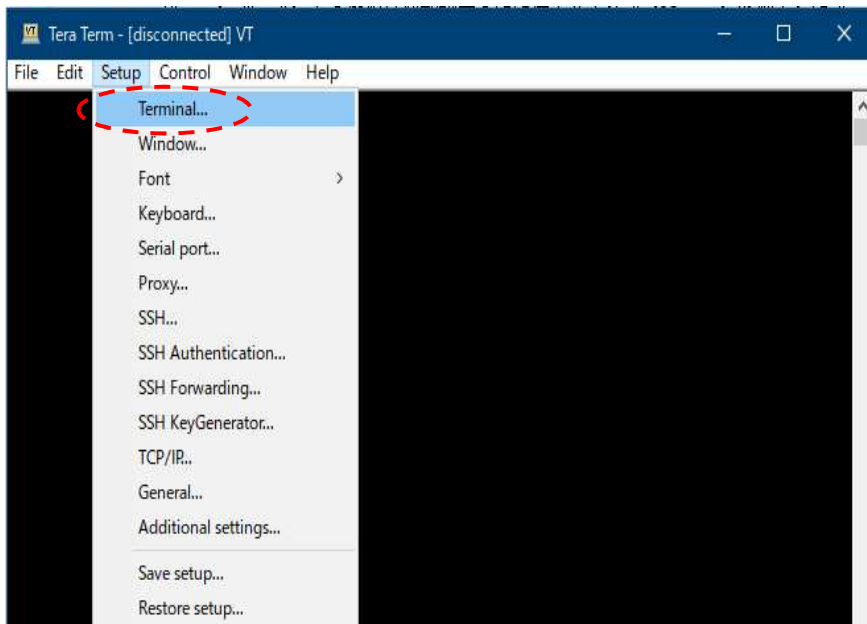



14. Input "K" and press the Enter key

Response data will be displayed as below if connection is OK.

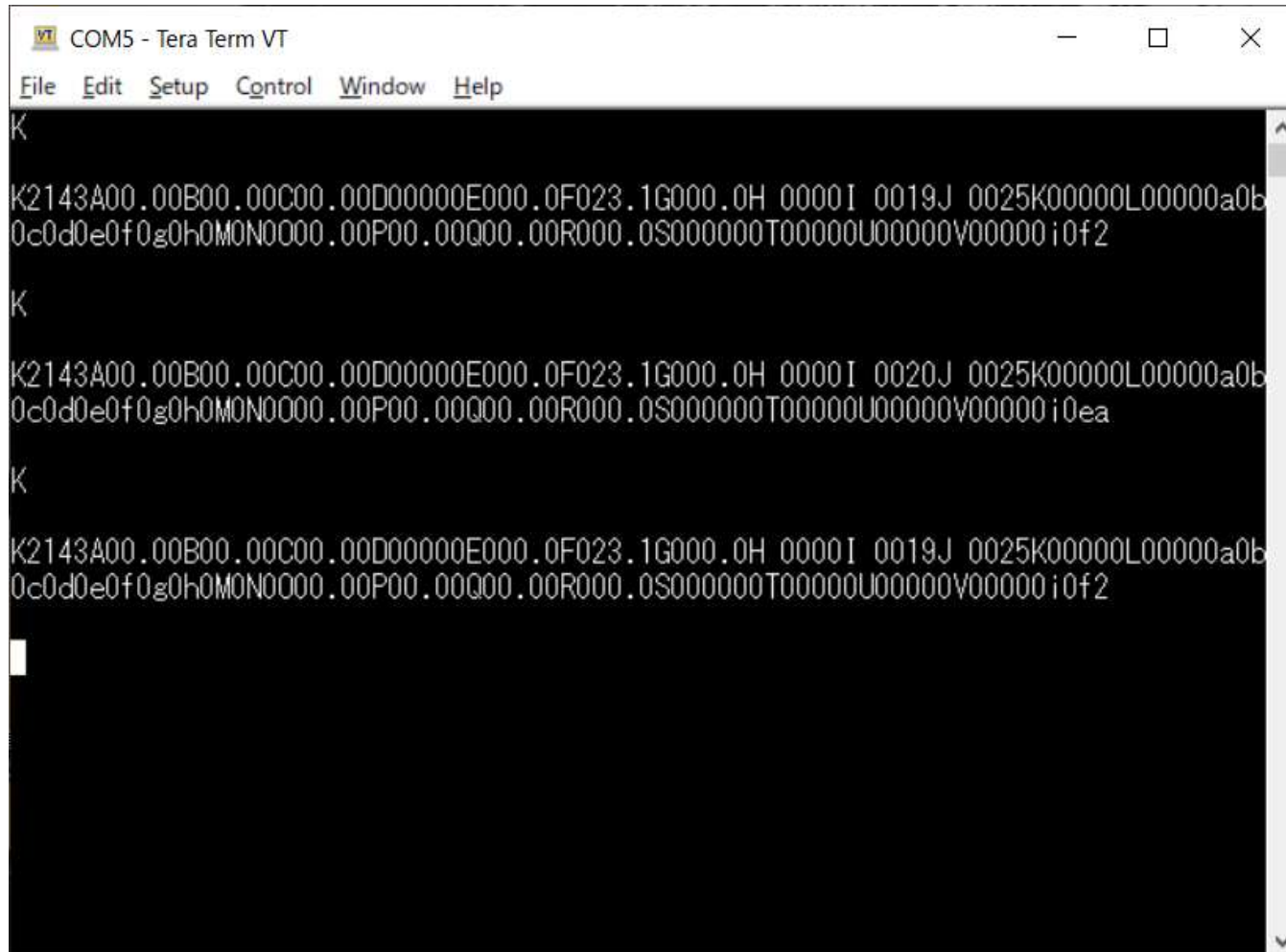

15. If the "KK" was displayed as the response from DBB-27, the machine did not output the data due to old SW version(before V3.0B), or improper setting/connection.